

How to change a sample with the Janis Cryofurnace and the Low Temperature Stick (5-500 K):

1. Remove old sample:
 - a. Disable DAS alarm by going to the Sample Environment computer on the KVM switch and pressing the “Mail Enabled” button on the Lakeshore_336 application (Figure 1).
 - b. Ensure that the sample mount and VTI temperature are between 100 and 300 K. (To speed cooling, once the furnace is below 325K, 200-300 mbar of He may be added.)
 - c. Turn off all heaters by hitting the red “ALL OFF” button on the front of the LakeShore 336 controller, (Figure 2). The adjacent red lights for Control Outputs 1 & 2 should turn off.
 - d. In the sample pit, unplug the 12-pin connector from the stick (Figure 3). The LakeShore alarm will beep.
 - e. Start flowing He, by turning the three-way valve to point to the gas line on the left (Figure 4).
 - f. Leaving He flowing, remove the three clamps.
 - g. Once the pressure on the gauge has reached approximately 900 mbar, carefully remove the sample stick, handling it by the black part, and hang it in the stick holder.
 - h. Cover the hole with the blank (Figure 5) and turn the three-way valve to vacuum on the right.
 - i. Once the pressure gauge is reading approximately 0 mbar, pump-purge the system by turning the three-way valve to He flow, then back to vacuum once the pressure gauge is reading approximately 900mbar.

Warning: Various sections of the sample stick may be hot or cold. Wear thermal gloves if needed to protect from extreme temperatures.
 - j. Remove the heat shield by removing the four screws. Be careful not to bend the stick.
 - k. Scan the sample with the RadEyeG radiation monitor, remove from stick and place with its barcode tag in the appropriate beamline location, on the wooden bench in the RMA.
2. Load new sample:
 - a. Use a heat gun and Kimwipes to completely dry any condensation on the stick. Do not leave the heat gun focused on one spot for a long time, to avoid damaging the wires.
 - b. Place the new sample onto the stick. Replace the heat shield, aligning the marks.
 - c. Turn the three-way valve to He.
 - d. Once the pressure on the gauge has reached approximately 900mbar, remove the blank, while continuing He flow.
 - e. Carefully place the sample stick into the cryofurnace, aligning the marks.
 - f. Turn the three-way valve to vacuum to the right.
 - g. Replace and tighten the three stick clamps.
 - h. Once the pressure gauge is reading approximately 0 mbar, turn the three-way valve to He flow.
 - i. Once the pressure gauge is reading approximately 900mbar, turn the three-way valve to vacuum.
 - j. Repeat steps j&k twice.
 - k. For high temperature (>300K) operation only, leave the sample space under vacuum.
 - l. For low-temperature (≤ 300 K) operation only, leave the sample space with 200-300 mbar of He to the cryofurnace by turning the three-way valve vertical to isolate the sample space once the desired pressure is reached.
 - m. Reattach the 12-pin plug to the sample stick.
 - n. On the front of the LakeShore 336 panel, reset the alarm by turning it off and immediately back on again. Consecutively press “Alarm +/-”, select *NO*, select *input C*, select *Alarm*, press up arrow to set it to *Off*, and press “Enter”. Then, select *Alarm* again, press down arrow to set it to *On*, and press “Enter”.
 - o. Turn the VTI heater on by pressing “A”, “Heater Range”, up arrow to select *High*, then “Enter”.
 - p. Turn the sample heater back on by repeating step q, starting with button “C” instead of “A”.
 - q. Re-enable DAS alarm by going to the Sample Environment computer on the KVM switch and pressing the “Mail Disabled” button on the Lakeshore_336 application.



Figure 1. Lakeshore_336 application on the Sample Environment computer

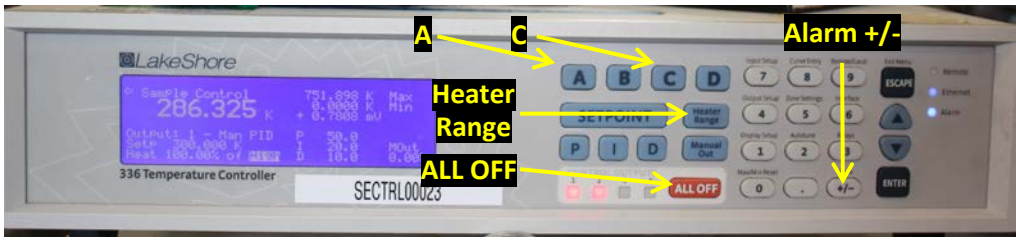


Figure 2. LakeShore336 controller

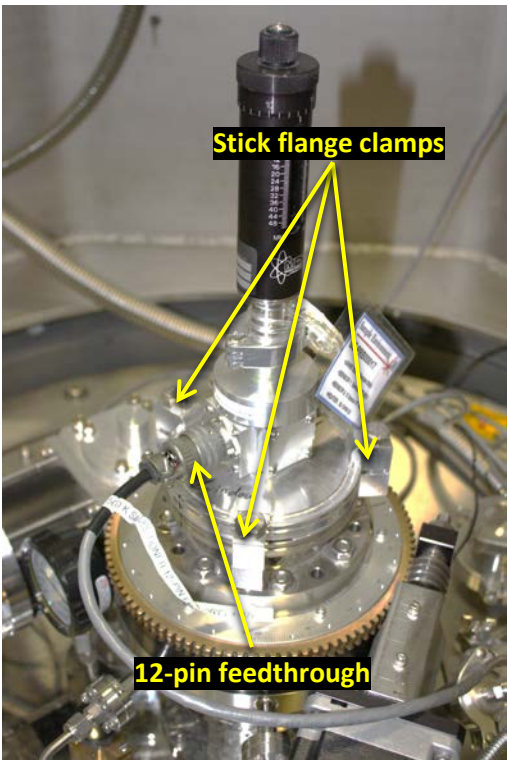


Figure 4. Sample stick connections

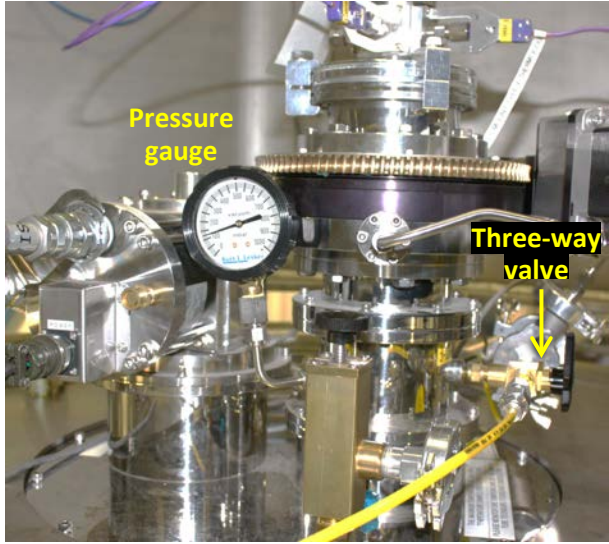


Figure 3. Pressure gauge and three-way valve



Figure 5. Flange blank in place on the cryofurnace