


# INSTRUMENT

BEAM LINE

# 1B

SPALLATION NEUTRON SOURCE


# Fact Sheet

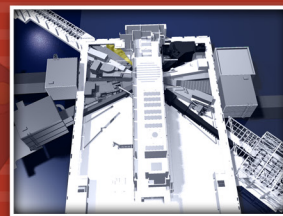

## NOMAD – NANOSCALE-ORDERED MATERIALS DIFFRACTOMETER

NOMAD is a high-flux, medium-resolution diffractometer that uses a large bandwidth of neutron energies and extensive detector coverage to carry out structural determinations of local order in crystalline and amorphous materials. The instrument enables studies of a large variety of samples, ranging from liquids and solutions, glasses, and nanocrystalline materials to long-range-ordered crystals. The enhanced neutron flux at SNS, coupled

with the advanced neutron optics and detector features, allows for unprecedented access to high-resolution pair distribution functions, small-contrast isotope substitution experiments, small sample sizes, and parametric studies.

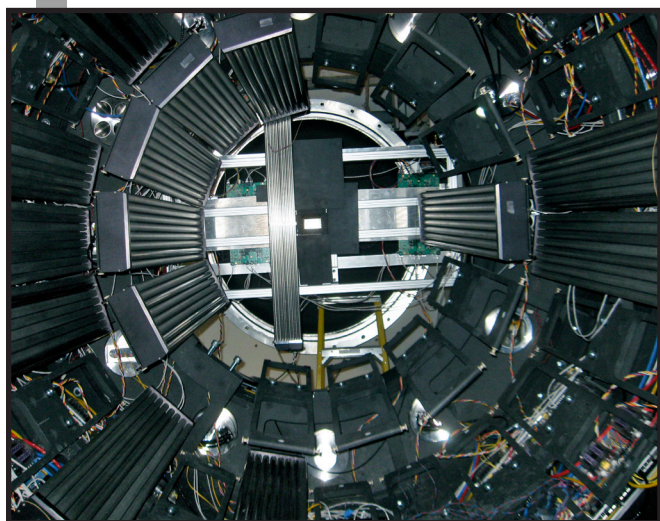

### SPECIFICATIONS

Moderator	Decoupled poisoned supercritical hydrogen
Moderator-to-sample distance	19.5 m
Sample-to-detector distance	0.5–3 m
Wavelength range	0.1–3 Å
Detector angular range	3–175° scattering angle
Initial coverage	4.0 sr
Full detector complement	8.2 sr
Flux on sample	$\sim 1 \times 10^8$ neutrons cm <sup>-2</sup> sec <sup>-1</sup>

Status: Available to users

### APPLICATIONS

- Environmental (e.g., solvent) effects on and direction of nanoscale structure formation
- In situ structural changes in nanoscale oxide catalysts used in automobile catalytic converters
- Structure of hydrogen storage materials under in situ conditions
- Transient structures of materials under extreme conditions (e.g., at high temperature or high pressure under the influence of transient fields or in metastable states)

### FOR MORE INFORMATION, CONTACT

Instrument Scientist: Jörg Neuefeind, [neuefeindjc@ornl.gov](mailto:neuefeindjc@ornl.gov), 865.241.1635

Instrument Scientist: Mikhail Feygenson, [feygensonm@ornl.gov](mailto:feygensonm@ornl.gov), 865.661.4813

[neutrons.ornl.gov/nomad](http://neutrons.ornl.gov/nomad)